


The Model **HLC-807** is a high performance liquid chromatograph designed specifically for preparative purpose. The system incorporates the many unique features developed by Toyo Soda that open up new possibilities for LC analysis in a number of fields. The adjacent chromatogram is an example of the 807's excellent performance.

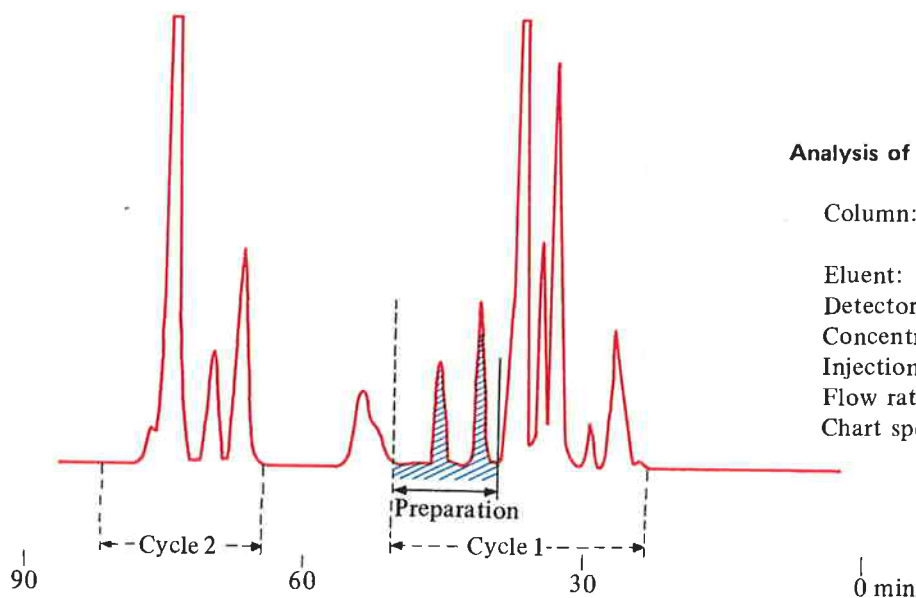

**Analysis of stabilizer for polyvinyl chloride**

Column: G 2000 HL  
(1 inch x 4 feet)  
Eluent: chloroform  
Detector: RI  
Concentration: 5%  
Injection volume: 100 $\mu$ l  
Flow rate: 6 ml/min  
Chart speed: 18 cm/hr